


# Slit Lamp Microscope

## S350C

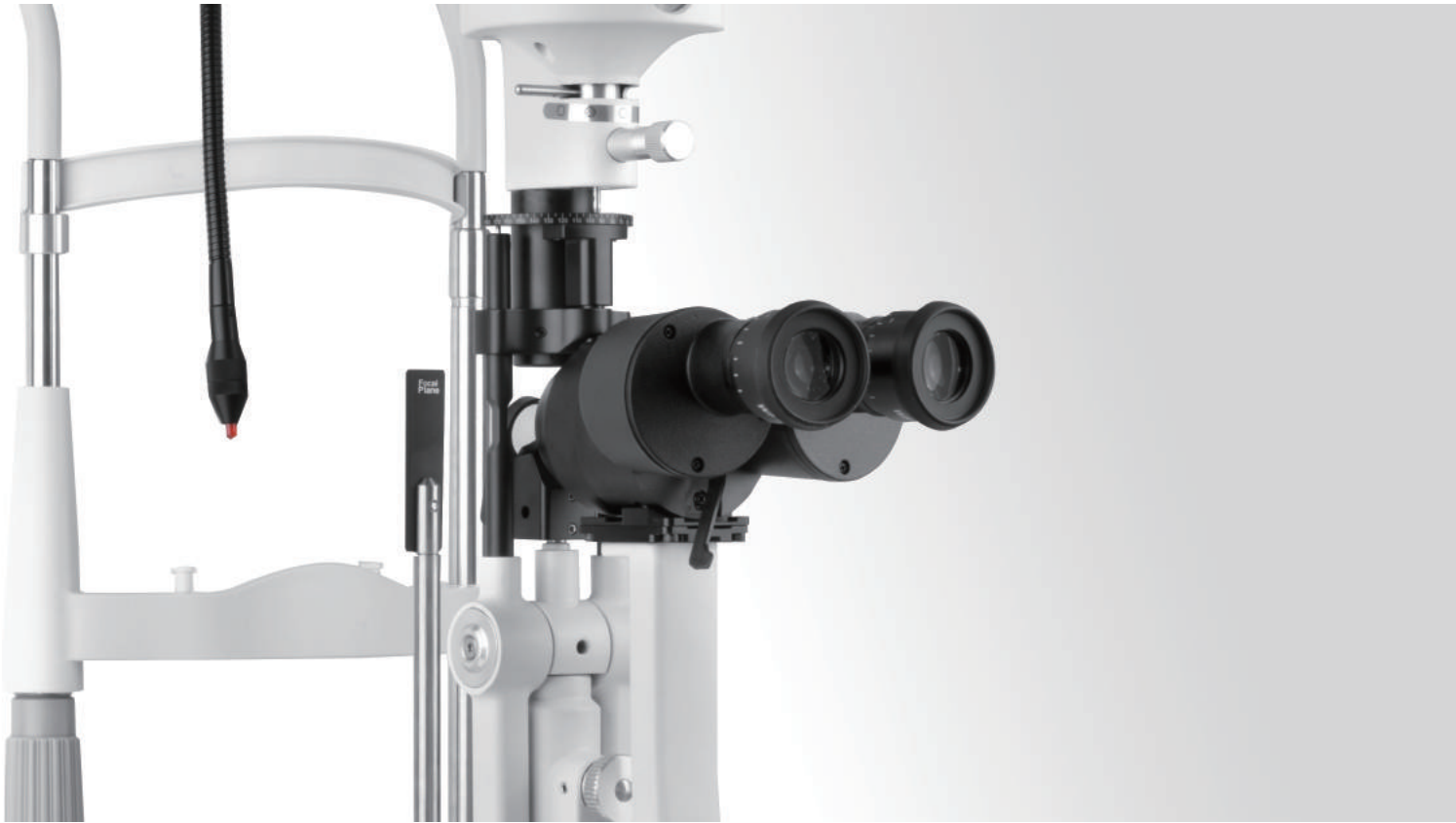

MediWorks adopts converging stereoscopic microscope used in the S350C. The excellent image for users. Two steps of magnification are provide.

The 14mm aperture fit for eye pre-checking of wearing contact lenses and the corneal foreign object examination.

6V/20W Osram halogen lamp is employed in S350C, which provides strong and stable lighting for observations. LED light source is optional.

Heat absorbing filter,ND filter,Red-free filter and Cobalt Blue filter are provided.

### Microscope

Microscope Type	Converging Stereoscope
Magnification Change	2 steps by Objective Lens Rotation
Total Magnification	10X, 16X
Eyepieces	10X
Angle between Eyepieces	13°
Pupillary Adjustment	52mm-78mm
Diopter Adjustment	+5D~-8D
Field of View	Ø18mm(10X), Ø14.5mm(16X)

### Slit Illumination

Slit Width	0~14mm continuous (slit becomes a circle at 14mm)
Slit Length	1~14mm continuous
Aperture Diameters	Ø14mm, Ø10mm, Ø5mm, Ø3mm, Ø2mm, Ø1mm, Ø0.2mm
Slit Angle	0°~180°
Slit Inclination	5°, 10°, 15°, 20°
Filters	Heat-absorbing filter, ND filter, Red-free filter, Cobalt Blue filter
Lamp	6V/20W Halogen Lamp
Luminance	≥50KLX

### Instrument Base

Longitudinal Movement	110mm
Lateral Movement	110mm
Vertical Movement	30mm
Fine Base Movement	15mm

### Chin-rest

Vertical Movement	80mm
Fixation	LED

### Power Supply

Input Voltage	220V/110V~±10%
Input Frequency	50Hz/60Hz
Power Consumption	30VA (max)
Output Voltage	Lamp 6V, Fixation 3V

### Packaging

Dimension	770mm x 470mm x 570mm(L/W/H)
Gross weight	24kgs
Net weight	16kgs



\*Subject to change in design or specification without advance notification (\*Table is optional)

Follow Us on



www.mediworks.com.cn

**Shanghai MediWorks Precision Instruments Co.,Ltd.**

Add:Zone A,2nd Floor,No.69,Lane 1985,Chunshen Road,Minhang District,  
Shanghai,200237,China

Tel: +86-21-54260421 54260423 Fax:+86-21-54260425

Email: international@mediworks.biz

Distributed by